Impella 5.5[®] with SmartAssist[®]

Circulatory Support System

USER MANUAL

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MANUALE D'USO

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KÄYTTÖOHJE

MANUAL DO UTILIZADOR

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ABIOMED[®]

Recovering hearts. Saving lives.

TABLE OF CONTENTS

English	1
Deutsch	17
Italiano	35
Français	53
Čeština	71
Español.....	87
Svenska.....	105
Dansk	123
Nederlands	139
Norsk.....	157
Polski.....	173
Suomi	191
Português.....	209
Magyar	227
Slovak.....	245

INTENDED AUDIENCE

This Instructions for Use and Clinical Reference Manual is designed for healthcare professionals. It contains clinical and technical information to guide healthcare professionals in their use of the Impella 5.5® with SmartAssist®. To use the system you must understand and follow these instructions. The Impella 5.5 with SmartAssist may be used only for its intended purpose.

WARNINGS

Warnings alert you to situations that can cause death or serious injury. The dark gray  symbol appears before warning messages.

Use of the Impella® System by trained and experienced practitioners has been associated with improved outcomes. Consequently, the first use of Impella should be preceded by the completion of a contemporary Abiomed Impella training program and include on-site proctoring during the first use by Abiomed clinical support personnel certified in the use of Impella.

Fluoroscopy is required for the insertion of the Impella guidewire and Impella 5.5 with SmartAssist Catheter.

During insertion, avoid manual compression of the inlet, outlet, or sensor areas of the cannula assembly on the Impella 5.5 Catheter.

The sterile components of the Impella 5.5 System can be used only if the packaging is not damaged, and the expiration date has not elapsed.

Do **NOT** resterilize or reuse Impella 5.5 with SmartAssist Catheter. It is a disposable device and is intended for single use only. Reuse, reprocessing, or resterilization may compromise the structural integrity of the catheter and/or lead to catheter failure which, in turn, may result in patient injury, illness, or death.

Retrograde flow will occur across the aortic valve if the Impella 5.5 with SmartAssist Catheter is set at P-level P-0.

Do **NOT** use saline in the purge system.

Do **NOT** use alcohol-based fluids for lubrication.

Do **NOT** use an Impella 5.5 with SmartAssist System if any part of the system is damaged.

To prevent the risk of explosion, do **NOT** operate the Impella 5.5 with SmartAssist System near flammable anesthetics.

If at any time during the course of support with the Impella 5.5 with SmartAssist Catheter, the Automated Impella Controller alarms "Purge Pressure Low" or "Purge System Open," follow the instructions presented in this manual.

 MR Unsafe - Do **NOT** subject a patient who has been implanted with an Impella 5.5 with SmartAssist Catheter to magnetic resonance imaging (MRI). The strong magnetic energy produced by an MRI machine may cause the Impella 5.5 System components to stop working, and result in injuries to the patient. An MRI may also damage the electronics of the Impella 5.5 System.

During defibrillation, do **NOT** touch the Impella 5.5 with SmartAssist Catheter, cables, or Automated Impella Controller.

Avoid overinserting the Impella 5.5 Catheter and possibly impinging the catheter tip against the walls of the vasculature or ventricles.

Do **NOT** advance or withdraw the Impella 5.5 with SmartAssist Catheter against resistance without using fluoroscopy to determine the cause of the resistance. Doing so could result in separation of the catheter or guidewire tip, damage to the catheter or vessel, or ventricular perforation.

In patients with transcatheter aortic valves position the Impella system carefully to avoid interaction with the transcatheter aortic valve prosthesis. Unintentional interaction of the Impella motor housing with the TAVR device may result in destruction of the impeller blades. This can lead to systemic embolization, serious injury, or death. In this situation, avoid repositioning while the device is running; turn the device to P0 during repositioning or any movement that could bring the outlet windows into proximity to the valve stent structures. If there is low flow observed in a patient implanted with a transcatheter aortic valve prosthesis, consider damage of the impeller and replace the Impella as soon as possible.

To reduce the risk of cardiac or vascular injury (including perforation) when manipulating the heart during cardiac surgery, evaluate the position of the pump using imaging guidance prior to manipulating the heart, and monitor position. In instances where the Impella pump has been placed prior to performing cardiac surgery with aortic cross clamping and cardioplegic arrest, care should be taken when manipulating the heart when the pump position is fixed with application of the aortic cross clamp across the pump catheter.

To reduce the risk of cardiac or vascular injury (including ventricular perforation) when advancing or torquing the Impella, adjustments should be performed under imaging guidance.

To reduce the possibility of fibers being drawn into the Impella, customers should avoid exposing the inlet and cannula section of the Impella Heart Pumps to any surfaces or fluid baths where the device can come into contact with loose or floating fibers.

To reduce the risk of cardiac injury (including ventricular perforation), physicians should exercise special care when inserting the Impella Catheter in patients with complex anatomy. This includes patients with known or suspected: decreased ventricular cavity size, ventricular aneurysms, congenital heart disease, or compromised cardiac tissue quality in the settings of acute infarction with tissue necrosis.

To reduce the risk of vascular injury, physicians should exercise caution when inserting the Impella Catheter in patients with complex peripheral vascular anatomy. This includes patients with known or suspected: unrepai red abdominal aortic aneurysm, significant descending thoracic aortic aneurysm, dissection of the ascending/ transverse/descending aorta, chronic anatomical changes in the relationship of the aorta/aortic valve/ventricular alignment, significant mobile atheromatous disease in the thoracic or abdominal aorta or peripheral vessels.

Physicians should exercise special care when inserting the Impella Catheter during active Cardiopulmonary Resuscitation (CPR). In addition, active CPR maneuvers may change the position of the Impella Device, introducing the risk of cardiac or vascular injury (including ventricular perforation). Check that the pump is positioned correctly after CPR with echocardiography guidance.

CAUTIONS

Cautions indicate situations in which equipment may malfunction, be damaged, or cease to operate. The light gray ⚠ symbol appears before caution messages.

Handle with care. The Impella catheter can be damaged during removal from packaging, preparation, insertion, repositioning and removal. Do **NOT** bend, excessively torque, pull, or place excess pressure on the catheter or mechanical components at any time.

Patients with aortic stenosis or other abnormal aortic valve performance may be compromised by the use of the Impella 5.5 with SmartAssist Catheter. Patients with aortic valve disease should be observed for aortic insufficiency.

Use only original accessories and replacement parts supplied by Abiomed.

To prevent device failure, do **NOT** start the Impella 5.5 with SmartAssist Catheter until the placement guidewire has been removed.

Do **NOT** remove the Impella 5.5 with SmartAssist Catheter over the length of the placement guidewire.

When replacing the purge cassette, the replacement process must be completed within 90 seconds. The Impella 5.5 with SmartAssist Catheter may be damaged if replacement takes longer than 90 seconds.

Do **NOT** kink or clamp any part of the Impella 5.5 with SmartAssist Catheter.

Have a backup Automated Impella Controller, purge cassette, and Impella 5.5 with SmartAssist Catheter available in the unlikely event of a device failure.

Use with Shockwave Coronary Intravascular Lithotripsy (IVL) Catheter at a distance of less than 20 mm between the optical sensor and IVL device may interfere with or damage the Impella optical sensor. Prior to IVL therapy pulsing, physicians should assess and verify this distance. If the placement signal is not displayed, monitor patient hemodynamics and confirm Impella position with imaging and motor current pulsatility. Loss of the placement signal does not impact Impella hemodynamic support.

POSSIBLE COMPLICATIONS (EU)

There are risks of complications with every procedure using a blood pump. These include among others:

- Acute renal dysfunction
- Aortic valve injury
- Cardiogenic shock
- Cerebral Vascular Accident/Stroke
- Myocardial infarction
- Renal failure
- Thrombocytopenia
- Hemolysis
- Bleeding
- Limb Ischemia
- Immune reaction
- Embolism, thrombosis
- Cardiac or vascular injury (including ventricular perforation)
- Positioning problems causing haemolysis or reduced haemodynamic support
- Infection and septicemia
- Dislocation of the pump
- Cardiovalvular injuries due to extreme movement of the suction cannula in relation to the cardiac valve or as a result of attachment by suction of the pump to the valve system following incorrect positioning
- Endocardiac injuries as a result of attachment of the pump due to suction
- Pump failure, loss of pump components following a defect
- Patient dependency on the pump after use for support

OVERVIEW

The Impella 5.5 Catheter is an intravascular microaxial blood pump that supports a patient's hemodynamic system. The Impella 5.5 Catheter is inserted via axillary artery cut-down through the artery and into the left ventricle. The Impella 5.5 Catheter may also be surgically inserted when there is access to the ascending aorta through a sternotomy or thoracotomy.

When properly positioned, the Impella 5.5 Catheter delivers blood from the inlet area, which sits inside the left ventricle, through the cannula, to the outlet opening in the ascending aorta. Physicians and device operators monitor the correct positioning and functioning of the Impella 5.5 Catheter on the display screen of the Automated Impella Controller.

This section describes the components of the Impella® Catheter and the Automated Impella Controller, as well as the accessory components.

REUSABLE SYSTEM COMPONENTS

The Impella 5.5 System consists of the following reusable components:

- Automated Impella Controller—provides the user interface, alarm indications, and portable battery
- Automated Impella Controller cart—for easy transport of the Automated Impella Controller

SINGLE-USE SYSTEM COMPONENTS

The Impella 5.5 with SmartAssist System also includes the following single-use components:

- Impella 5.5 with SmartAssist Catheter
- Purge cassette
- 0.018 inch, 260 cm placement guidewire
- Impella Axillary Insertion kit
- Silicone Plugs
- Incision template

INTENDED USE (EU)

INTENDED USE IN THE EUROPEAN UNION

The Impella 5.5® with SmartAssist® heart pump is an intracardiac pump for supporting the left ventricle. It is intended for clinical use in cardiology and in cardiac surgery for up to 30 days for the following indications, as well as others:

- The Impella 5.5 with SmartAssist is a cardiovascular support system for patients with reduced left ventricular function, e.g., post-cardiotomy, low output syndrome, cardiogenic shock after acute myocardial infarction.
- The Impella 5.5 with SmartAssist may also be used as a cardiovascular support system during coronary bypass surgery on the beating heart, particularly in patients with limited preoperative ejection fraction with a high risk of postoperative low output syndrome.


CONTRAINDICATIONS (EU)

CONTRAINDICATIONS IN THE EUROPEAN UNION

The Impella 5.5® with SmartAssist® heart pump is contraindicated for the following situations:


- Mechanical aortic valves, severe aortic valvular stenosis
- Hematological disorder causing fragility of the blood cells or hemolysis
- Hypertrophic obstructive cardiomyopathy (HOCM)
- Aneurysm or necrotomy or severe anomaly of the ascending aorta and/or the aortic arch
- Mural thrombus in the left ventricle
- Ventricular septal defect (VSD) after myocardial infarction
- Anatomic conditions precluding insertion of the pump

SYSTEM CONFIGURATION

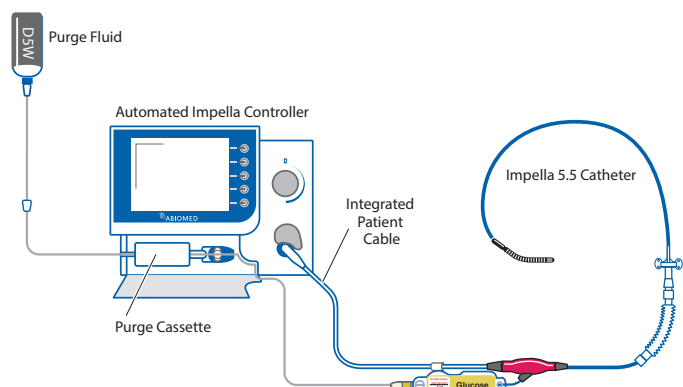IMPELLA 5.5® WITH SMARTASSIST® CATHETER

The Impella 5.5 with SmartAssist Catheter is an intravascular microaxial blood pump that delivers up to 5.5 liters of blood per minute from the left ventricle into the aorta.

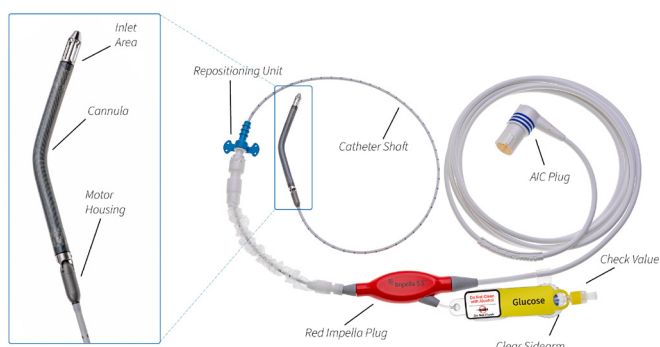PURGE CASSETTE

Do not use saline in the purge system.

The purge cassette delivers rinsing fluid to the Impella 5.5® with SmartAssist® Catheter. The purge fluid (typically 5% glucose solution in water with heparin) flows from the purge cassette through the catheter to the microaxial blood pump to prevent blood from entering the motor.


ACCESSORIES


The 0.018 inch, 260 cm guidewire is used for the placement of the Impella 5.5 Catheter. The guidewire has a radiopaque, shapable tip.

These alternative guidewires have been tested and approved for use with the Impella System: Boston Scientific Platinum Plus™ ST 0.018 in, and Boston Scientific V-18 Control Wire™ ST 0.018 in.

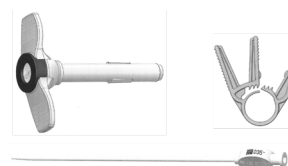The Impella 5.5 Catheter Set contains a 23 Fr diameter x 6 cm length peel-away introducer and two (2) graft locks used to attach a graft onto the introducer. (Note: Only one graft lock is required when used with the recommended Hemashield™ Platinum graft; a back-up is provided.)

The kit is packaged with an 8 Fr silicone-coated lubrication dilator. It is recommended to be used in conjunction with a 10 mm diameter x 20 cm length Hemashield Platinum graft.

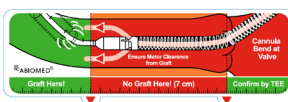The incision template is used to assist in correct placement of the Impella 5.5 Catheter during direct aortic insertion by indicating the correct location for the incision in the ascending aorta.

The two silicone plugs can be placed around the catheter shaft to help control bleeding during and after Impella 5.5 Catheter insertion.

Hospital Provided:

Solution (typically 5% glucose in water with 25 or 50 IU/mL of heparin) is used as the purge fluid through the Impella 5.5 Catheter.

The Automated Impella Controller cart holds the Automated Impella Controller. The cart has wheels for easy transport of the controller and a storage basket. (For more information, including assembly instructions, refer to the Automated Impella Controller cart instructions for use.)

PRE-SUPPORT EVALUATION

Before initiating the procedure, evaluate the patient for factors that may prevent successful placement of the Impella Catheter. Use imaging technology to examine the patient's vasculature access site. An echo assessment of the left ventricle is also recommended to rule out left ventricular thrombus, mechanical aortic valves, or severe aortic insufficiency. Consider inserting a diagnostic Pulmonary Artery (PA) catheter to provide continuous hemodynamic monitoring, including pulmonary artery and central venous pressures, measurement of cardiac output, and SvO₂.

STARTUP

SUPPLIES NEEDED

- Automated Impella Controller
- Impella 5.5® with SmartAssist® Catheter and accessories
- Diagnostic catheter (AL1 or MP without side holes or pigtail without side holes)
- 500 cc bag of glucose solution for purge solution in water (5% recommended; 5% to 20% acceptable) with 25 or 50 IU heparin/mL
- 10 mm x 20 cm Hemashield Platinum vascular graft (if using Axillary Insertion kit)
- Impella Axillary Insertion kit for axillary insertion of the Impella

Note: If the purge solution is supplied in bottles, open the vent on the purge fluid spike and follow the same procedure as if supplied in bags.

CASE START


To avoid fibers drawn into the Impella:

- Keep the Impella Heart Pump in its packaging tray until just before insertion.
- Do not attempt to run the pump in a basin of saline prior to insertion.
- Do not attempt to rinse and reinsert the device after initial insertion.
- Hold the surgical towel or 4 x 4 gauze pad away from the inflow and outflow windows, when controlling blood splatter during insertion of the Impella Heart Pump through the introducer.

If you miss a step in the process of setting up the Impella Catheter, or if you exceed the amount of time allowed to complete a step, the Automated Impella Controller will display additional instruction screens for continuing the setup process for some steps. For the steps discussed below, all shaded steps require sterile technique.

1. Press the **START NEW CASE** soft button from the startup screen or plug in a new Impella Catheter. "Case Start" can also be selected by pressing the **MENU** soft key.

2. The controller displays the screen shown below.

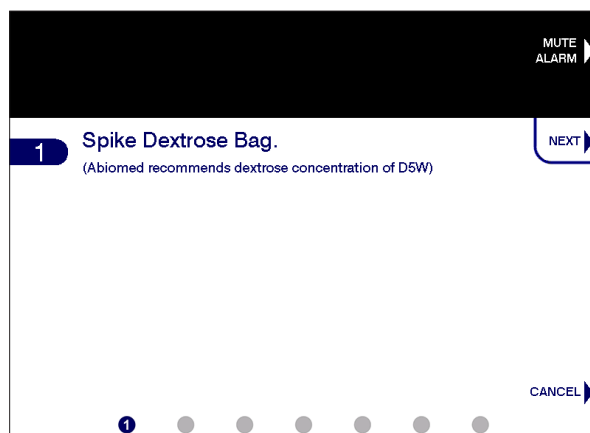
INSERT PURGE CASSETTE

1. Open the purge cassette package.
2. Secure the yellow luer connector on the purge tubing to the sterile field.
3. Pass the purge cassette and spike off the sterile field.
4. Spike the glucose bag/bottle.
5. Press the **NEXT** soft button to continue.
6. Open the purge cassette door by pressing the release on the left side of the controller. Insert the purge cassette into the Automated Impella Controller.
7. Insert the purge cassette into the compartment on the front of the controller. Follow the diagram on the inside of the purge cassette door for proper placement.
8. Slide the purge disc into the slot to the right of the purge cassette until it snaps into place. The controller will automatically begin priming the purge cassette.
9. Extend the purge tubing and close the purge cassette door. There is sufficient room around the edges of the purge cassette door so that it will not pinch the purge tubing as it exits.
10. The controller automatically begins priming the purge cassette after it is inserted. The progress bar marks the progress of the purge cassette priming.


CONNECT THE IMPELLA CATHETER

1. Remove the Impella Catheter from its package using sterile technique and inspect the catheter for damage.
2. Inspect the cable for damage, including damage to the connector pins at the controller end.
3. Pass the sterile connector cable from the Impella Catheter off the sterile field.


4. Open the cover on the blue catheter plug by rotating clockwise. Line up the notch on the connector cable with the notch in the blue catheter plug on the front of the Automated Impella Controller and plug the cable into the controller.

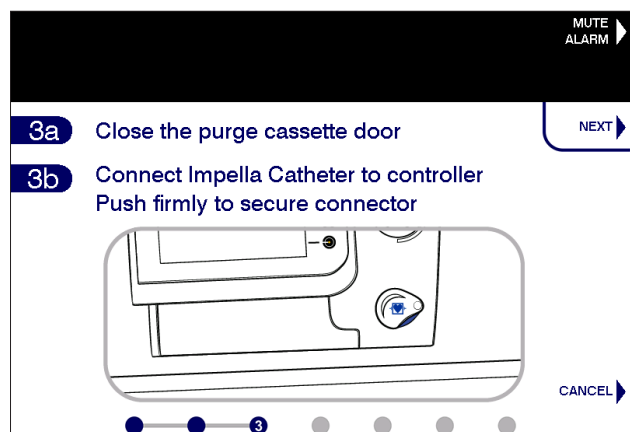
5. Once the purge cassette is primed and the controller detects that the connector cable is plugged in, it prompts you to connect the luer to the Impella Catheter.

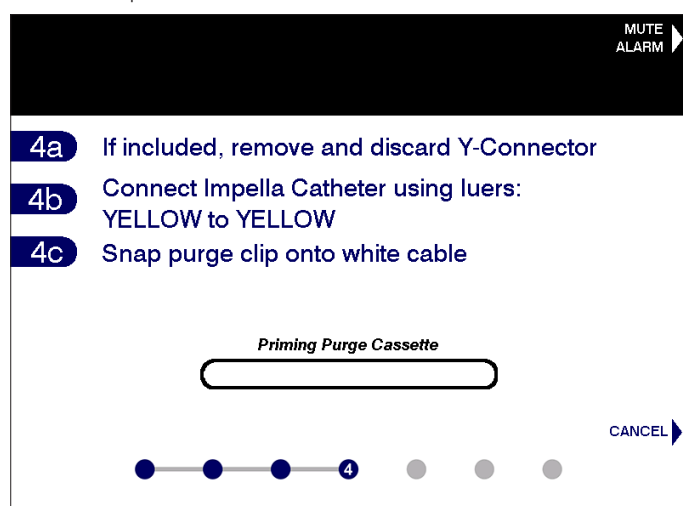6. Connect and tighten the yellow luer on the purge tubing to the Impella Catheter sidearm if not already complete.

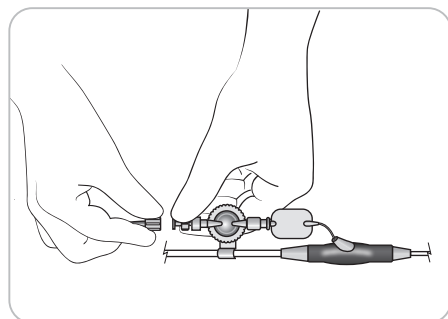
7. When the controller detects that the luer is connected, it automatically begins priming the purge lumen.

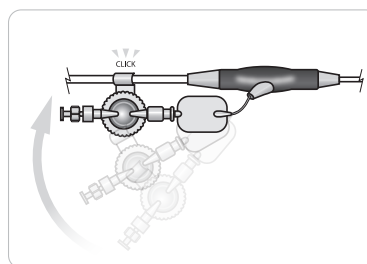8. Snap the purge clip (located on the pressure reservoir of the clear sidearm) to the connector cable
9. Once the purge lumen is primed, the controller automatically advances to the next screen.
10. The first step on the next screen prompts you to enter the purge fluid information.

ENTER PURGE FLUID DATA

1. Enter the purge fluid information.

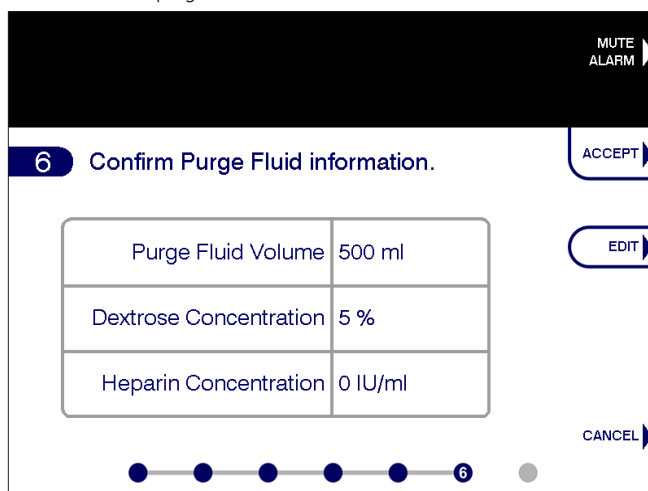2. To select the default values displayed on the screen, press the **ACCEPT** soft button. This will select those values and automatically advance to the next screen. Note: The default purge fluid values will be the purge fluid values from the last Case
3. To change the purge fluid information, press the **EDIT** soft buttons, scroll to the appropriate item and push the selector knob to select it or use the white arrow soft buttons. Then scroll through the values and push the selector knob or press **SELECT** to make a new selection. Press the **DONE** button to finish editing. The controller will use the default values if no other selections are made.
 - Purge fluid can be set to 50 mL, 100 mL, 250 mL, 500 mL, or 1000 mL.
 - Glucose (Dextrose) concentration can be set to 5%, 10% or 20%.
 - Heparin concentration can be set to 0 IU/ml, 5 IU/ml, 6.25 IU/ml, 10 IU/ml, 12.5 IU/ml, 20 IU/ml, 25 IU/ml, 40 IU/ml, 50 IU/ml.

SECURE THE PURGE TUBING

1. To complete the setup, connect the purge tubing to the white connector cable by pushing the purge tubing into the clips attached to the white connector cable.

AXILLARY INSERTION OF THE IMPELLA 5.5® CATHETER

NOTE – Proper surgical procedures and techniques are the responsibility of the medical professional. The described procedure is furnished for information purposes only. Each physician must evaluate the appropriateness of the procedure based on his or her medical training and experience, the type of procedure, and the type of systems used.

If the Impella Catheter is used in the OR as part of open heart surgery, manipulation may be performed only at the access site. Direct manipulation of the catheter assembly through the aorta or ventricle may result in serious damage to the Impella Catheter and serious injury to the patient.

When using a pigtail diagnostic catheter with side holes, ensure that the guidewire exits the end of the catheter and not the side hole. To do so, magnify the area one to two times as the guidewire begins to exit the pigtail.


During placement of the Impella Catheter, take care to avoid damage to the inlet area while holding the catheter and loading the placement guidewire.

-  The introducer and graft lock are supplied sterile and can be used only if the packaging is not damaged and the expiration date has not elapsed.
-  Fluoroscopy is required for the insertion of the Impella guidewire and Impella Catheter.
-  During insertion, avoid manual compression of the inlet or outlet areas of the Impella Catheter or the sensor area of the cannula on the Impella 5.5 with SmartAssist Catheter.
-  The graft must be affixed to the introducer proximal to the retainers on the introducer repositioning unit to prevent the introducer from sliding out of the graft.
-  When inserting the Impella Catheter through the introducer and into the graft, be sure to clamp the graft with a vascular clamp just above the anastomosis to avoid blood loss through the pump cannula during insertion through the valve.
-  The Impella Axillary Insertion kit is intended to be used for insertion only. To provide continued hemostasis, the introducer must be peeled away and the repositioning unit inserted into the graft.
-  Do **NOT** resterilize or reuse any components of the Impella Axillary Insertion kit. All components are disposable and intended for single use only. Reuse, reprocessing, or resterilization may compromise performance.


-  The introducer is designed to be inserted into a graft. It is not intended for direct insertion into the artery.
-  Abiomed recommends the use of a 10 mm diameter Hemashield Platinum graft with the introducer for proper fit and hemostasis between the graft and the introducer. A smaller diameter graft may not fit over the introducer.
-  Abiomed recommends the use of a 20 cm length graft to allow enough length to fully insert the Impella Catheter cannula into the graft prior to releasing vascular clamps at the anastomosis to minimize blood loss through the cannula.
-  Do **NOT** kink or clamp the white portion of the Impella Catheter with cm markings with anything other than a soft jaw vascular clamp. Do **NOT** kink or clamp the peel-away introducer.
-  Proper positioning of the Impella Catheter is extremely important and it is worthwhile to take extra time when positioning the catheter.
-  Take care to insert the guidewire with diagnostic catheter into the middle of the hemostatic valve of the introducer to avoid tearing the valve.
-  When inserting the Impella Catheter into the introducer, take care to insert it straight into the center of the introducer valve.
-  Be sure that the stopcock on the peel-away introducer is always kept in the closed position. Significant bleed back can result if the stopcock is open.

The following steps describe the recommended technique for axillary artery insertion of the Impella 5.5® with SmartAssist® Catheter.

1. Isolate and expose the axillary artery and obtain control via proximal and distal vessel loops.
2. Attach a 10 mm diameter x 20 cm long vascular graft to the axillary artery using a standard end-to-side anastomosis. **NOTE:** Abiomed recommends using a Hemashield Platinum graft and recommends using at least a 60 degree bevel on the end of the graft to facilitate passage of the rigid motor housing into the artery.
3. Clamp the graft with a vascular clamp just above the anastomosis and loosen the vessel loops to allow blood to flow into the graft to assess for hemostasis at the anastomosis.
4. Insert the introducer into the graft and secure it with one (1) provided graft lock. To place the graft lock, open it and place it between the retainers and the hub on the introducer to prevent the introducer from sliding out of the graft. **NOTE:** If a graft other than the Hemashield Platinum is used, two graft locks may be required to maintain hemostasis between the graft and the introducer. Correct positioning of the second graft lock is illustrated in below.

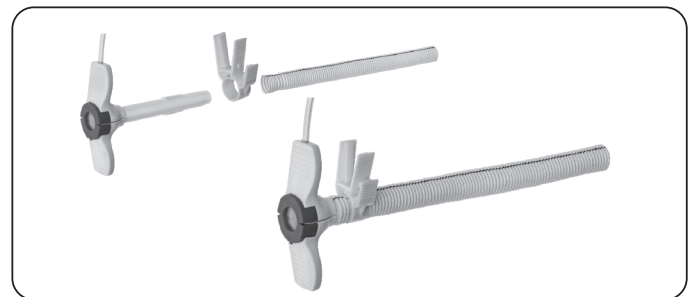
Introducer, Graft Lock, and Hemashield Platinum Graft (Graft Not Supplied)

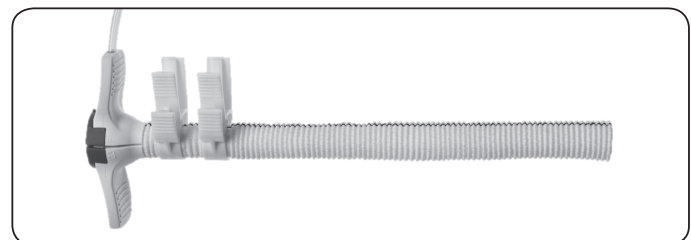Correct Positioning If Second Graft Lock Required


5. Secure the graft lock by pressing both the outside tabs together. When fully closed, the graft lock provides hemostasis. If hemostasis is not achieved, make sure to press the two tabs together to fully close the graft lock as shown below. The graft lock cannot be damaged by over closing. **NOTE:** The graft may also be secured over the introducer using heavy sutures or umbilical tape.

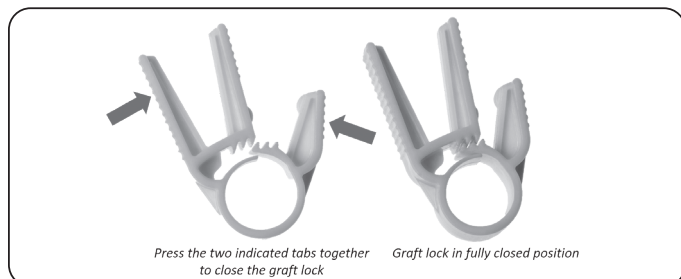Closing the Graft Lock


6. Remove the vascular clamp on the graft and insert a 0.035 inch diagnostic guidewire with a 4–6 Fr diagnostic catheter into the introducer, taking care to center the wire and catheter in the center of the hemostatic valve. Advance the guidewire and catheter into the left ventricle.
7. Remove the diagnostic guidewire and exchange it for a stiff 0.018 inch placement guidewire. With the 0.018 inch placement guidewire properly positioned in the left ventricle, remove the diagnostic catheter.
8. Administer heparin per institutional protocol. When the ACT is greater than or equal to 250 seconds, remove the dilator.
9. Remove the protective sleeve on the provided 8 Fr silicone-coated lubrication dilator, being careful to avoid getting silicone on your hands. Insert the dilator into the introducer over the 0.018 inch placement guidewire to coat the hemostatic valve with silicone oil to facilitate insertion of the Impella Catheter through the hemostatic valve assembly. Once fully inserted, remove the dilator, keeping the 0.018 inch placement guidewire in place.
10. Clamp the graft with a vascular clamp just above the anastomosis to avoid blood loss through the pump cannula during insertion through the valve.
11. While maintaining guidewire position, backload the Impella Catheter onto the 0.018 inch placement guidewire and advance the catheter over the guidewire through the introducer into the graft such that the entire pump cannula and motor housing resides in the graft and only the catheter shaft is seen exiting the valve.
12. Remove the vascular clamp and continue inserting the Impella Catheter into the aorta. Continue advancing across the aortic valve using fluoroscopic imaging to properly position the cannula bend at the aortic valve annulus, placing inlet approximately 5 cm deep into ventricle. Remove the placement guidewire and initiate Impella Catheter support as described later in this section.
13. Clamp the graft adjacent to the axillary artery with a soft jawed vascular clamp or have an assistant apply digital pressure to control bleeding at the base of the graft so that the introducer can be removed and the graft shortened. **NOTE:** To ensure the soft jaw vascular clamp is completely sealing over the graft and the 9 Fr catheter, open the sidearm flush valve on the introducer and verify blood is not leaking from the system.
14. Slide the repositioning sheath back to the red Impella plug. **NOTE:** For the 23F x 6cm peel-away introducer, it may be necessary to pull the introducer over the repositioning sheath to remove the peel-away completely from the artery.

15. To remove the introducer, release the graft lock by pressing the two adjacent long tabs together as shown below and remove it from the graft.

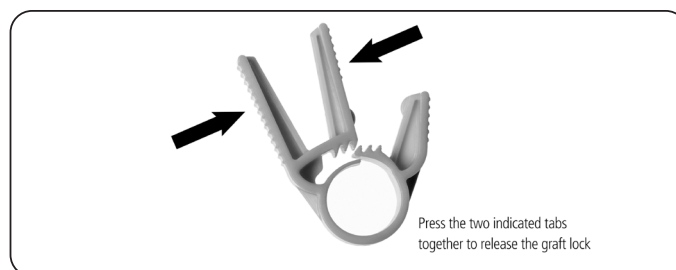Releasing the Graft Lock


16. Slide the introducer fully out of the graft prior to peeling it away. To peel the introducer off the catheter shaft, crack the hub by applying pressure to the thumb tabs and then peel the sheath off the catheter. **NOTE:** When breaking the hemostatic valve in the sheath hub, the valve may stretch before separating.
- a). Grasp the two “wings” and bend back until the valve assembly comes apart. Continue to peel the two wings until the introducer is completely separated from the catheter shaft.

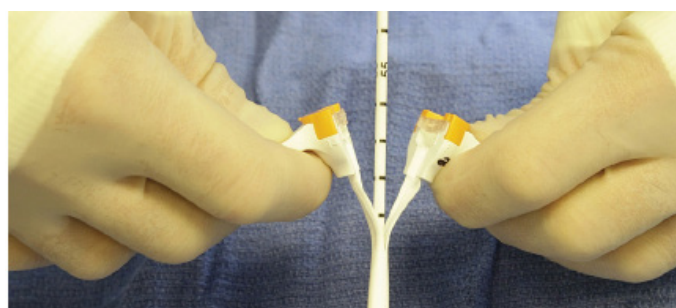17. Trim any excess graft and slide the blue suture hub into the graft. **NOTE:** The hub should be at the skin level and the length of the remaining graft material should be just long enough to secure the graft around the blue suture hub with all of the graft buried beneath the skin.
18. Using heavy silk suture, secure the graft around the blue suture hub so that the position of the Impella Catheter can still be adjusted. Remove the vascular clamp adjacent to the axillary artery.
19. The wound should be closed over the trimmed graft with the end of the blue suture hub clearly visible. Anchor the hub securely to the skin.
20. If there is slack in the catheter, remove the excess slack. Verify placement with fluoroscopy and with the placement signal. If the Impella 5.5 Catheter advances too far into the left ventricle and the controller displays a ventricular waveform rather than an aortic waveform, follow steps a-c below.
 - a) Pull the catheter back until an aortic waveform is present on the placement screen.
 - b) When the aortic waveform is present, pull the catheter back an additional 3cm for Impella 5.5 with SmartAssist. (The distance between adjacent markings on the catheter is 1 cm.)
 - c) The catheter should now be positioned correctly.
21. Attach the anticontamination sleeve to the blue suture hub. Lock the anchoring ring in place by turning it clockwise. Secure the catheter shaft in place by tightening the Tuohy-Borst. Extend the sterile sleeve to maximum length and secure the end closest to the red Impella® plug by tightening the anchoring ring.

ALTERNATE INSERTION TECHNIQUE USING A SIDARM GRAFT & SILICONE PLUGS

1. After exposing the axillary artery and making the incision as described in the steps above, prepare a Dacron® vascular graft (10 mm x 20 cm) by beveling the end of the graft at a 45 to 60 degree angle.
2. Tighten the distal and proximal vessel loops to control bleeding.
3. Attach the vascular graft using the standard end-to-side anastomosis

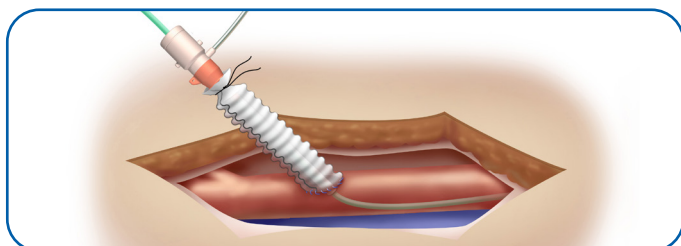4. Assess the anastomosis for hemostasis.
5. Attach a standard 6 Fr introducer to the distal end of the graft.
6. Advance a diagnostic catheter (Abiomed recommends a 6 Fr AL1 or Multipurpose without side holes or 4–5 Fr pigtail with or without side holes) over a diagnostic 0.035 inch or 0.038 inch guidewire into the left ventricle.
7. Remove the diagnostic guidewire and exchange it for the supplied 0.018 inch placement guidewire.
8. Tighten the vessel loops to control bleeding and remove the 6 Fr introducer.
9. Moisten the Impella 5.5 Catheter and push one of the silicone plugs onto the catheter shaft adjacent to the Impella 5.5 Catheter motor.
10. Backload the Impella 5.5 Catheter onto the 0.018 inch guidewire (as described in the steps in previous step 11).
11. With the graft held at the base, place the Impella 5.5 Catheter into the open end of the graft up to the level of the silicone plug.
12. Secure umbilical tape around the silicone plug.
13. Loosen both vessel loops and advance the Impella 5.5 Catheter along the guidewire into the left ventricle until it is properly positioned.

To prevent device failure, do not start the Impella 5.5 Catheter until the placement guidewire has been removed.

Do NOT remove the Impella 5.5 Catheter over the length of the placement guidewire.

14. Remove the guidewire.
15. Place a soft-jawed clamp at the anastomosis and remove the silicone plug.
16. Trim any excess graft and slide the repositioning sheath into position.
17. Using a heavy silk tie or umbilical tape, secure the graft around the blue hub of the repositioning sheath.
18. Close the wound over the trimmed graft with the end of the blue suture hub clearly visible. The steering catheter for the Impella 5.5 can be manipulated if needed by unsecuring the repositioning unit and moving the catheter in or out.
19. If there is slack in the catheter, remove the excess slack. Verify placement with fluoroscopy and with the placement signal. If the Impella 5.5 Catheter advances too far into the left ventricle and the controller displays a ventricular waveform rather than an aortic waveform, follow steps a-c below.

a) Pull the catheter back until an aortic waveform is present on the placement screen.

b) When the aortic waveform is present, pull the catheter back an additional 3cm for Impella 5.5 with SmartAssist. (The distance between adjacent markings on the catheter is 1 cm.)

c) The catheter should now be positioned correctly.

20. Attach the anticontamination sleeve to the blue suture hub. Lock the anchoring ring in place by turning it clockwise. Secure the catheter shaft in place by tightening the Tuohy-Borst. Extend the sterile sleeve to maximum length and secure the end closest to the red Impella® plug by tightening the anchoring ring.

DIRECT AORTIC INSERTION

NOTE – Proper surgical procedures and techniques are the responsibility of the medical professional. The described procedure is furnished for information purposes only. Each physician must evaluate the appropriateness of the procedure based on his or her medical training and experience, the type of procedure, and the type of systems used.

Avoid manual compression of the inlet, outlet, or sensor areas of the cannula assembly.

Do **NOT** kink or clamp the Impella Catheter with anything other than a soft jaw vascular clamp. Do **NOT** kink or clamp the peel-away introducer.

Handle with care. The Impella catheter can be damaged during removal from packaging, preparation, insertion, repositioning and removal. Do **NOT** bend, excessively torque, pull, or place excess pressure on the catheter or mechanical components at any time.

An incision larger than 6 mm may allow the front plug to advance into the aorta.

The Impella 5.5 Catheter is surgically inserted when there is access to the ascending aorta through a sternotomy or thoracotomy. Transesophageal echocardiography (TEE) is required to guide placement.

INSERTION PREPARATION

When using a pigtail diagnostic catheter with side holes, ensure that the guidewire exits the end of the catheter and not the side hole. To do so, magnify the area one to two times as the guidewire begins to exit the pigtail.

1. Using the supplied sterile incision template for positioning, place a sidebiter clamp on the aorta at least 7 cm above the valve plane. An incision too close to the aortic valve annulus could result in the catheter outlet area in the graft rather than the aorta.
2. Make an incision (or punch) no larger than 6 mm at the insertion site on the ascending aorta. The incision must be ≤ 6 mm in length to prevent the front silicone plug from advancing into the aorta through the incision.
3. Attach the Dacron® vascular graft (10 mm x 15 cm) to the aorta using the standard end-to-side anastomosis.
4. Administer heparin and achieve ACT of at least 250 seconds. **Note:** Maintaining ACT at or above 250 seconds will help prevent a thrombus from entering the catheter and causing a sudden stop on startup. *If the patient is receiving a GP IIb/IIIa inhibitor, the Impella 5.5 Catheter can be inserted when ACT is 200 or above.
5. When the anastomosis is complete, place a clamp at the distal end of the graft and then release the proximal clamp at the base of the graft. Examine the suture line for leaks and reclamp the graft at the base.

- Moisten the Impella 5.5 Catheter and push both silicone plugs up against the motor housing as shown below.

Impella 5.5 Catheter with Silicone Plugs

- With the graft clamped at the base, place the Impella 5.5 Catheter into the open end of the graft up to the level of the rear plug.
- When the catheter is in position, secure a tourniquet around the rear silicone plug. Tighten the tourniquet sufficiently to control bleeding around the rear plug while still allowing the catheter to slide through the plug.
- Release the clamp and advance the Impella 5.5 Catheter into the aorta.
- If the patient is on cardiopulmonary bypass (CPB), allow the heart to fill by restricting the return flow to the bypass machine and reducing CPB flow to a minimum setting, as long as acceptable physiologic systemic flow is maintained.
- As soon as the motor housing has passed into the aorta, use a ligature to loosely secure the front silicone plug flush to the graft. The silicone plug should be in the most proximal portion of the graft. **Note:** There should be no movement of the front silicone plug within the graft; however, the catheter shaft should move without resistance within the plug. When securing the front silicone plug to the graft, do not penetrate the silicone plug too deeply as this could cause damage to the Impella 5.5 Catheter.
- While the catheter is being advanced in the aorta, the initial placement signal has the characteristics of an aortic placement signal. Do not allow the front plug to advance beyond the base of the graft.
- To aid in passing the catheter through the aortic valve, apply slight pressure to the posterior aspect of the aortic valve to produce temporary aortic insufficiency.
- Gently advance the catheter forward until the inlet crosses the aortic valve and the bend of the catheter is at the level of the aortic valve annulus. Confirm with TEE guidance.

POSITIONING AND STARTING THE IMPELLA 5.5® WITH SMARTASSIST® CATHETER


When the Impella 5.5 with SmartAssist Catheter is not correctly placed, there is no effective unloading of the ventricle (hydraulic short circuit). The patient may not be benefiting from the flow rate shown on the controller.

Retrograde flow will occur across the aortic valve if the flow rate of the Impella Catheter is less than 0.5 L/min.

- Reconfirm that the placement guidewire has been removed and Confirm that the controller displays a pulsatile waveform and the cannula bend at the aortic valve annulus, placing inlet approximately 5 cm deep into ventricle. (See step 6 if the controller displays a ventricular waveform.)

- Press **START IMPELLA** soft button to open the P-level menu. Turn the selector knob to increase the P-level from P-0 to P-2.

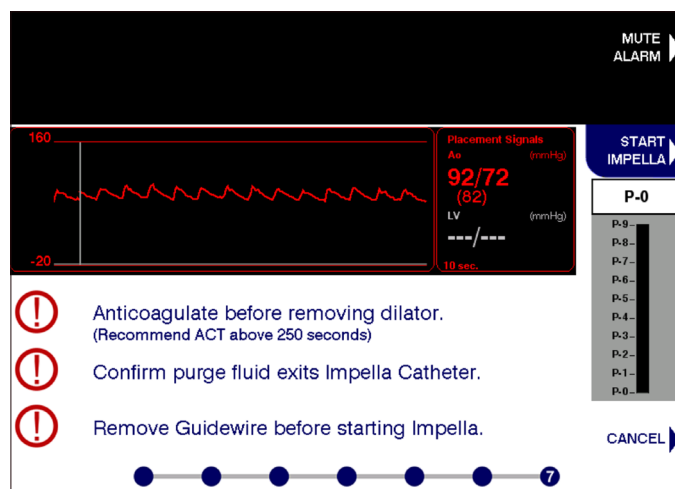- Press the selector knob to select the new P-level. Increase the P-level to P-9 to confirm correct and stable placement.
- Evaluate the catheter position in the aortic arch. If there is slack in the catheter, loosen the Touhy-Borst, remove the excess slack, and RE-TIGHTEN THE TUOHY-BORST. The catheter should align against the lesser curvature of the aorta rather than the greater curvature. Verify placement with fluoroscopy and with the placement signal screen.
- Reposition the catheter as necessary.
- If the Impella Catheter advances too far into the left ventricle and the controller displays a ventricular waveform (shown above) rather than an aortic waveform, follow these steps to reposition the catheter.
 - Loosen the Touhy-Borst, pull the catheter back until an aortic waveform is present on the placement screen.
 - When the aortic waveform is present, pull the catheter back an additional 3 cm for Impella 5.5 with SmartAssist. (The distance between adjacent markings on the catheter is 1 cm.)
 - RE-TIGHTEN THE TUOHY BORST
 - The catheter should now be positioned correctly.

USE OF THE REPOSITIONING UNIT

- Slide the blue suture hub and advance the ribs into the graft.
- Secure the repositioning unit to the patient with the blue suture pads or a StatLock® stabilization device.
- Evaluate the catheter position in the aortic arch. If there is excess slack, loosen the Touhy-Borst, remove the excess slack, and RE-TIGHTEN THE TUOHY-BORST. The catheter should align against the lesser curvature of the aorta rather than the greater curvature. Verify placement with fluoroscopy and with the placement signal.
- Attach the anticontamination sleeve to the blue section of the repositioning unit. Lock the anchoring ring in place by turning it clockwise. Secure the catheter shaft in place by tightening the connected anchoring ring.
- Carefully extend the anticontamination sleeve to maximum length and secure the end closest to the red Impella plug by tightening the anchoring ring.
- Select the lowest P-level that will enable you to achieve the highest flow rate necessary for patient support. You can select one of ten P-levels (P-0 to P-9) for the Impella 5.5 Catheter.

P-LEVEL

In **P-LEVEL** mode you can select one of ten P-levels (P-0 to P-9) for the Impella 5.5 with SmartAssist Catheter (see table below). Select the lowest P-level (P-2 or higher) that will enable you to achieve the highest flow rate necessary for patient support.

P-level	Mean Flow (L/min) 30 - 60 mmHg	Revolutions Per Minute (rpm)
P-0	0	0
P-1	0	12,000
P-2	0.0 - 1.9	17,000
P-3	1.1 - 2.7	20,000
P-4	1.9 - 3.3	22,000
P-5	2.8 - 3.7	24,000
P-6	3.4 - 4.1	26,000
P-7	3.9 - 4.5	28,000
P-8	4.3 - 4.9	30,000
P-9	5.0 - 5.5	33,000

**Flow rate can vary due to suction or incorrect positioning*

To operate the Impella Catheter in P-level mode:

1. Press the **FLOW CONTROL** soft button to open the **FLOW CONTROL** menu.
2. Turn the selector knob to increase or decrease the flow rate.
3. Press the selector knob to select the new flow rate.

ADJUSTING THE LV PLACEMENT SIGNAL

The LV placement signal and LV estimates are not displayed when pumps are running at P-3 or lower. Increase pump speed to P-4 or higher to re-enable signal

Alarm conditions and low pump speeds may affect the LV placement signal and estimates

Disruption of the outlet pressure placement signal, including alarms related to the Ao placement signal, will prevent calculation and display of an LV estimate. An operational Ao placement signal is required for the LV estimate


LV placement signal calibration is not available if the P-level is less than P-4 or if Suction, Placement Signal Not Reliable, or positioning alarms are active

Abnormal conditions, including cardiac arrhythmia, Ao-LV uncoupling, or aortic stenosis may limit the utility of the LV adjustment tool

ADJUST LV PLACEMENT SIGNAL


Adjust the LV placement signal to reduce potential measurement variability. Adjust the LV Placement Signal when the white notification appears on the screen. The notification appears first when a suggested adjustment is calculated. If a suggested adjustment is calculated, then a second notification will appear after 24 hours of pump use.

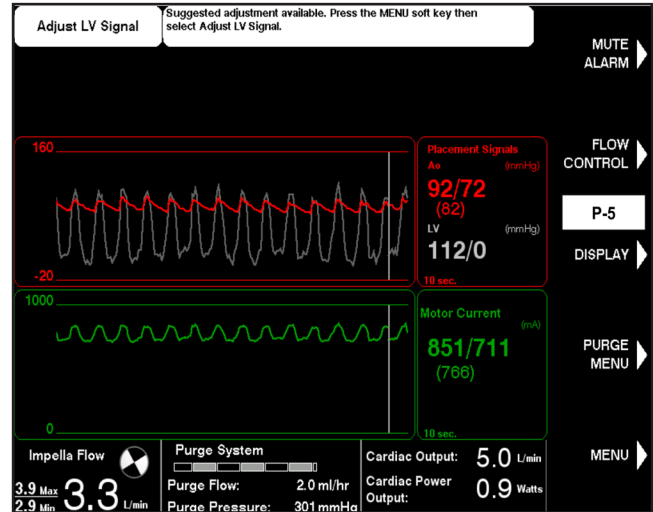
To adjust LV Placement Signal:

1. Press the **MENU** soft buttons.
2. Select "Adjust LV Signal" option with rotary knob
3. Adjust the waveform using toggle arrows or rotary knob. The LV Adjustment defaults to the suggested adjustment value.
4. Press **DONE** to confirm suggested adjustment.
LV waveform adjustments occur in increments of 1 mmHg from -60 to 60 mmHg.

Recommended: Do not adjust to a value other than the suggested value.

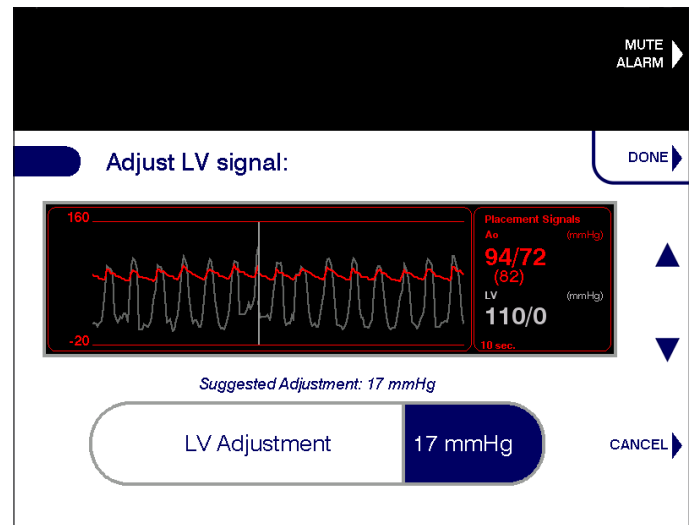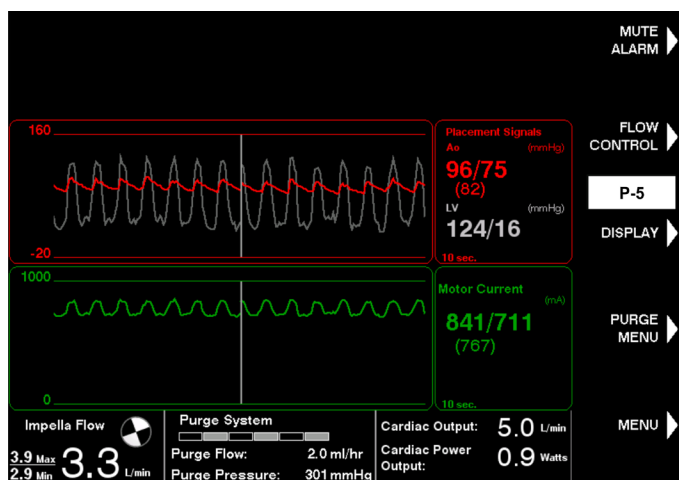
Note: It is not atypical to see a LV systolic value higher than the Ao systolic

ENTERING CARDIAC OUTPUT


The cardiac output and cardiac power output metrics are for informational purposes only. Do not use for diagnostic purposes or patient monitoring. Independently verify all parameters displayed with a cleared or approved diagnostic device. Enter cardiac output into the Automated Impella Controller from a reference device, such as a Swan-Ganz catheter. Enter a new cardiac output every 8 hours. After 7 hours from entry, a white notification will trigger to enter a new cardiac output.

TO ENTER CARDIAC OUTPUT

1. Press the **MENU** soft button.
2. Select the "Enter Cardiac Output" using the rotary knob.
3. Enter the total Cardiac Output, which can be any value from 0.0 to 10.0 L/min in increments of 0.1 L/min
4. Press the **DONE** soft button or the rotary knob to complete.

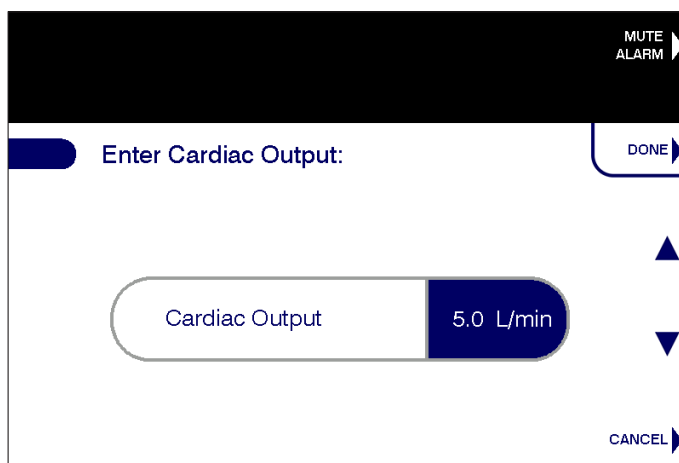A Cardiac Output Confirmation will be displayed if a total Cardiac Output is entered that is less than or equal to the current Impella Flow. The **CONFIRM** soft button will use the current Impella Flow as the total Cardiac Output entry and the Native Cardiac Output will not be trended. The **BACK** soft button will go back to the Enter Cardiac Output screen. The **CANCEL** soft button will exit the workflow.

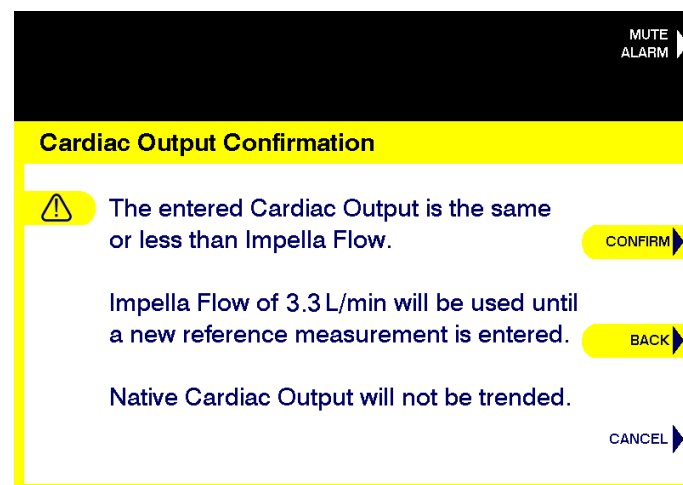CARDIAC OUTPUT, NATIVE CARDIAC OUTPUT & CARDIAC POWER OUTPUT CALCULATIONS

Once a cardiac output is entered, the Automated Impella Controller can calculate an initial cardiac power output and native cardiac output using the following equations:

$$CPO = (CO \times MAP) / 451$$

$$NCO = CO - \text{Impella Flow}$$

The Native Cardiac Output estimate is derived from a relationship between native function and aortic pulse pressure (PP). This relationship is linear and scaled by a calibration factor, β which may vary between patients and as an individual patient's condition changes. This relationship can be shown by the following equation:

$$NCO = \beta \times PP$$

Once the calibration factor is obtained, the Automated Impella Controller will continue to calculate the cardiac output, native cardiac output, and cardiac power output for the next 8 hours using the above equations.

If cardiac power output values are below or equal to 0.6, the value will display as yellow.

Note: Do not use values as a clinical diagnostic tool, this is for informational purposes only.

CARDIAC OUTPUT ENTRY REMINDERS

A white reminder to enter the CO will appear every 15 minutes for the last hour prior to the 8 hour timeout. If a new CO is not entered after 8 hours, the values will show as dash marks until a new cardiac output is entered.

A white reminder will also be displayed to update the CO if the Automated Impella Controller detects a significant change in the vascular state. This notification will be triggered if the average NCO or PP diverges from their initial values by a significant amount.

PURGE CASSETTE PROCEDURES

When replacing the purge cassette or purge fluid, the replacement process should be completed within 90 seconds of disconnecting the luer. The Impella Catheter may be damaged if replacement takes longer than 90 seconds of disconnecting the luer.

There are three procedures for maintaining the Impella Catheter purge system:

- Change cassette and purge fluid bag
- Change purge fluid bag
- De-Air purge system

Each procedure can be accessed using the **PURGE MENU** soft button. The purge cassette procedures are discussed below.

CHANGE CASSETTE AND BAG

The Purge Cassette may need to be replaced for extended use of the Impella system.


Follow these steps to change both the purge cassette and purge fluid:

1. Press **PURGE MENU** and select "Change Cassette and Bag" from the menu.
2. Select **START** to begin the cassette and fluid change process.
3. When prompted by the controller, disconnect the luer from the Impella catheter.

Note on Connecting Purge Tubing to the Catheter:

If included, remove and discard the Y-connector and connect the yellow luer on the purge tubing directly to the yellow check valve on the Impella Catheter.

4. Open the purge cassette door by pressing the button on the left side of the console. Remove and discard the old cassette and purge fluid bag.
5. Open the new purge cassette. Spike the new purge fluid bag with the new purge cassette tubing. Select **NEXT** to continue.
6. Insert the new purge cassette into the controller. Be sure to slide the purge disc into place and extend the purge tubing through the gap in the purge cassette door when you close the door.

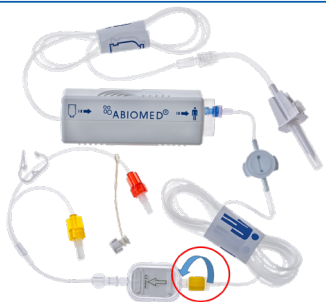7. Confirm the luer is disconnected. Press **NEXT** to proceed to prime the purge cassette.
8. Update the purge fluid information.
 - a) To select the default purge fluid values displayed on the screen, select **CONFIRM**.
 - b) To change the purge fluid information, select **EDIT**. Then use the soft keys to navigate selections and edit values. Select **DONE** to complete editing.
9. When Steps 1 through 8 are complete, connect the luer from the new purge cassette to the Impella catheter.

CHANGE PURGE FLUID BAG

Note on Purge Solution Bottles:

If the purge solution is supplied in bottles, open the vent on the purge fluid spike and follow the same procedure as if supplied in bags.

These are the steps you will follow to change only the purge fluid:

1. Press **PURGE MENU** and select "Change Purge Fluid Bag."
2. Select **START** to begin the purge fluid change process.
3. When prompted by the controller, remove the old purge bag and replace by spiking the new purge fluid bag. Select **NEXT** to advance to the next step.
4. Update the purge fluid information.
 - a) To select the default purge fluid values displayed on the screen, select **CONFIRM**.
 - b) To change the purge fluid information, select **EDIT**. Then use the soft keys to navigate selections and edit values. Select **DONE** to complete editing.
5. When prompted by the controller, disconnect the luer from the Impella Catheter. The controller will automatically prime the tubing, which will flush the fluid from the last bag out of the purge cassette tubing.
Note: the instructions to disconnect the luer and to automatically prime the tubing only occurs if the user changed the purge fluid concentration
 - a) To skip the flush select **SKIP PRIME**.
6. When prompted by the controller, connect the yellow luer from the purge cassette to the Impella catheter.


DE-AIR PURGE SYSTEM

These are the steps you will follow to de-air the purge system:

1. Press **PURGE MENU** and select "De-Air Purge System."
2. Select **START** to begin the de-air process.
3. Make sure that the purge fluid bag is NOT empty or inverted and that the tubing is NOT kinked. Select **NEXT** to continue.
4. Disconnect the purge tubing from the Impella Catheter.
5. Confirm that no air remains in the purge tubing. If air remains, press **BACK** to repeat the air removal process.
6. Connect the purge tubing to the luer on the Impella Catheter to complete the de-air procedure.

AIR DETECTED ALERT

During any of the purge system processes above, the controller automatically monitors for air in the system. If air is detected in the system, the controller provides an alert to disconnect the luer, as shown below. Once the luer is disconnected, the controller automatically de-airs the purge system.

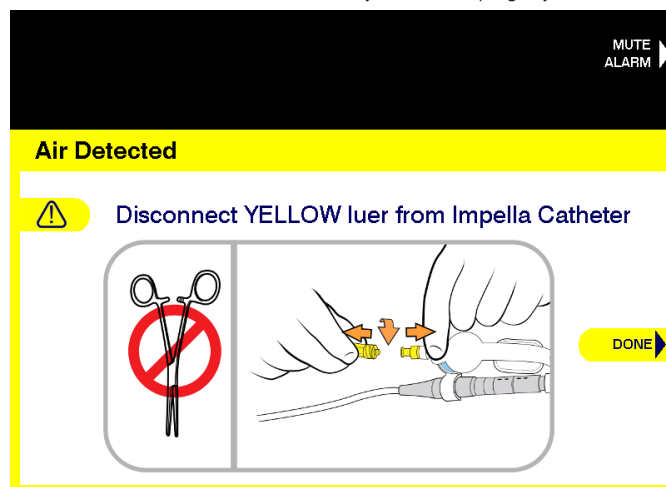PATIENT WEANING

Weaning the patient from the Impella Catheter is at the discretion of the physician. The Impella 5.5® with SmartAssist® has been approved for ≤ 30 days of use. However, weaning could be delayed beyond the normal use for temporary support as an unintended consequence of continued instability of the patient's hemodynamics. Inability to wean the patient from the device within a reasonable time frame should result in consideration of a more durable form of left ventricular support. The LVEDP/CO Trend can be used while weaning. Refer to the Automated Impella Controller Instructions for Use document for more information.

The following weaning instructions are provided as guidance only.

RAPID WEANING

1. Initiate rapid weaning by decreasing catheter P-level in 2-level steps at intervals of several minutes (for example, P-6 to P-4 to P-2). Do NOT decrease P-level to below P-2 until just before removing the catheter from the ventricle.
2. When the P-level has been reduced to P-2, maintain the patient on P-2 support for at least 10 minutes before discontinuing circulatory support.
3. If the patient's hemodynamics remain stable, decrease the P-level to P-1, pull the catheter into the aorta, and stop the motor by decreasing the P-level to P-0.
4. Explant the catheter.
5. Follow institutional guidelines for arterial closure.
6. Disconnect the connector cable from the Automated Impella Controller and turn the controller off by pressing the power switch on the side of the controller for 3 seconds.

SLOW WEANING

1. Initiate slow weaning by decreasing catheter P-level in 2-level steps over time as cardiac function allows (for example, P-6 to P-4 to P-2). Do NOT decrease P-level to below P-2 until just before removing the catheter from the ventricle.
2. When the P-level has been reduced to P-2, maintain the patient on P-2 support until the patient's hemodynamics remain stable before discontinuing circulatory support.
3. If the patient's hemodynamics remain stable, decrease the P-level to P-1, pull the catheter into the aorta, and stop the motor by decreasing the P-level to P-0.
4. Explant the catheter.
5. Follow institutional guidelines for arterial closure.
6. Disconnect the connector cable from the Automated Impella Controller and turn the controller off by pressing the power switch on the side of the controller for 3 seconds.

REMOVING THE IMPELLA 5.5® WITH SMARTASSIST® CATHETER

The Impella Catheter can be removed after weaning when the introducer is still in place or when the catheter is secured with the repositioning sheath.

REMOVING THE IMPELLA CATHETER WITH THE INTRODUCER IN PLACE

1. Reduce the P-level to P-0.
2. Remove the Impella Catheter through the introducer.
3. Wait until ACT drops below 150 seconds.
4. When ACT is below 150 seconds, remove the introducer.

5. Disconnect the connector cable from the Automated Impella Controller and turn the controller off by pressing the power switch on the side of the controller for 3 seconds.
6. Apply manual compression for 40 minutes or per hospital protocol.

REMOVING THE IMPELLA CATHETER SECURED WITH THE REPOSITIONING UNIT

1. When ACT is below 150 seconds, press **FLOW CONTROL** and reduce the P-level to P-0.
2. Remove the Impella Catheter and repositioning unit together (the catheter will not come through the repositioning unit).
3. Disconnect the connector cable from the Automated Impella Controller and turn the controller off by pressing the power switch on the side of the controller for 3 seconds.
4. Apply manual compression for 40 minutes or per hospital protocol.

VASCULAR GRAFT CLOSURE

When closing the vascular graft, consider individual patient characteristics and select the strategy most consistent with an optimal clinical result. The entire vascular graft can be removed if indicated, but it is not mandatory to do so. Graft closure options include:

- Amputating the vascular graft and sewing the small end-to-side remnant closed by hand
- Using a vascular stapler to close the graft near the surface of the aorta
- Removing the complete graft with local patch closure, if necessary

OPERATING THE IMPELLA 5.5 WITH SMARTASSIST CATHETER WITHOUT HEPARIN IN THE PURGE SOLUTION

The Impella Catheter is designed to be operated with a purge solution that contains heparin. Operation of the system without heparin in the purge solution has not been tested. In the event that a patient is intolerant to heparin, due to heparin-induced thrombocytopenia (HIT) or bleeding, physicians should use their clinical judgment to assess the risks versus benefits of operating the Impella System without heparin.

If it is in the best interest of the patient to operate the system without heparin, the dextrose solution is still required, and physicians should consider *systemic delivery* of an alternative anticoagulant. The Impella Catheter has not been tested with any alternative anticoagulants in the purge solution. Use of alternative anticoagulants may reduce the longevity or performance of the Impella catheter.

OPERATING THE IMPELLA 5.5 WITH SMARTASSIST CATHETER IN ELECTROMAGNETIC FIELDS

The Impella Catheter contains a permanent magnet motor that emits an electromagnetic field. This field may produce electromagnetic interference with other equipment. In addition, other equipment that emits a strong electromagnetic field may affect the operation of the Impella Catheter motor.

ELECTROANATOMIC MAPPING (EAM) SYSTEMS

The electromagnetic field emitted by the Impella Catheter may produce interference with the magnetic location detection component of the electroanatomic mapping (EAM) system, particularly when the mapping catheter is close to the Impella Catheter motor. For example, mapping in the right or left ventricular outflow tracts places the mapping catheter in close proximity to the Impella Catheter motor in the ascending aorta.

- Electromagnetic interference may appear as:
- Instability in the displayed location of the mapping catheter
 - Magnetic interference errors generated by the electroanatomic mapping system

When operating the Impella® Catheter in the presence of an EAM system, use P-level mode. Operate the Impella Catheter at P-1–P-5 or P-7. The motor speeds at these P-levels cause the least interference. Best performance is observed when the Impella Catheter motor is at least 3 cm from the sensors in the mapping catheter. If you suspect interference, follow the troubleshooting steps below.

Observation	Actions
Interference with the magnetic location detection component of the EAM system	1. Check for and address other sources of interference. 2. Reposition the Impella Catheter to ensure that the Impella motor is at least 3 cm from the sensors in the mapping catheter; however do NOT pull the inlet area out of the left ventricle. 3. Ensure that the Impella Catheter is operating at P-1–P-5 or P-7, as these P-levels cause the least interference.

MAGNETIC NAVIGATION SYSTEMS (MNS)

- When initiating Impella Catheter support in the presence of a magnetic navigation system (MNS), follow the steps below:
1. Insert the Impella Catheter following the steps outlined in this manual.
 2. Place the MNS magnets in the “Reduced” or “Stowed” position.
 3. Start the Impella Catheter in the manner described in this manual. Increase P-level to at least P-5.
 4. Place the MNS magnets in the “Navigate” position and proceed with magnetic navigation.

Keep operating the Impella Catheter at a P-level of at least P-5 when the MNS magnets are in the “Navigate” position. If the P-level falls below P-5, the Impella Catheter may stop running.

To resume operation, follow the steps below.

During magnetic navigation of the mapping catheter, the motor current of the Impella Catheter may temporarily increase to the point that the catheter stops running. The table below explains how to resume operation.

When the MNS magnets are in the “Navigate” position, the displayed Impella Catheter flow may be artificially elevated. To accurately assess the flow rate, note the displayed flow when the magnets are in the “Stowed” position.

Observation	Actions
Unable to start Impella or Impella stops running	1. Place the MNS magnets in the “Reduced” position and attempt to start the Impella Catheter. 2. If the Impella Catheter does NOT start with the magnets in the “Reduced” position, place the magnets in the “Stowed” position and start the Impella Catheter. 3. Increase the Impella Catheter P-level to P-5 or higher. 4. Place the MNS magnets in the “Navigate” position and proceed with magnetic navigation.
MNS magnets: “Navigate” Displayed flow seems too high or MNS magnets: “Stowed” Displayed flow drops	The Impella Catheter displayed flow will be artificially elevated when the MNS magnets are in the “Navigate” position. The displayed flow will be accurate when the MNS magnets are in the “Stowed” position.

ALARMS OVERVIEW

The Automated Impella Controller monitors various functions to determine whether specific operational parameters are within expected limits. When a parameter goes outside of its specified limits, the controller sounds an alarm tone and displays an alarm message that can be viewed on the display screen on the front of the controller. The alarm tone indicates the severity of the alarm. The alarm message on the display screen is color-coded for severity and provides details on the cause of the alarm and how to resolve the alarm.

For a full list of all Impella alarms, refer to the Automated Impella Controller Instructions for Use document.

ALARM LEVELS

Alarms are divided into three levels of severity:

- Advisory (white)
- Serious (yellow)
- Critical (red)

Category	Description	Audible Indicator*	Visual Indicator
Advisory	Notification	1 beep every 5 minutes	Alarm header on white background
Serious	Abnormal situation. Prompt action needed.	3 beeps every 15 seconds	Alarm header on yellow background
Critical	High priority. Immediate action needed.	10 beeps every 6.7 seconds	Alarm header on red background

* Sound pressure of audible alarm indicators is 60-80 dBA at a distance of 1 m.

MUTE ALARM FUNCTION

Pressing the **MUTE ALARM** button on the upper right of the Automated Impella Controller display screen will silence the audible alarm indicator for 2 minutes (for red or yellow alarms) or 5 minutes (for white advisory alarms). When an alarm is silenced, the words “MUTE ALARM” next to the button are replaced by the mute alarm indicator, a crossed-out bell icon.

The audible indicator will shut off if an alarm condition is resolved before you press **MUTE ALARM**. The visual message, however, will continue to be displayed, with the alarm header on a gray background, for 20 minutes or until you press **MUTE ALARM**. This allows you to identify the alarm that occurred.


SYMBOLS

	Caution; follow instructions for use
	Defibrillator-proof type CF equipment
	Keep dry
	Storage temperature (e.g. 10°C to 30°C)
	Declares conformity with Directive 93/42/EEC for medical devices or with Regulation (EU) 2017/745 on medical devices, and with Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment
2019-11-06	Date of manufacture (eg, 2019-11-06)
	Protect from sunlight
	Symbol for lot designation; the manufacturer's lot designation must be stated after the LOT symbol
123456	Abiomed part number (eg, part number 123456)
123456	Manufacturer's serial number (eg, serial number 123456)
2019-11-06	Use-by date (eg, use before 2019-11-06)
	Do not reuse
	Sterilized using ethylene oxide
	Medical Device
	Federal (USA) law restricts this device to sale by or on order of a physician: Prescription only.
	Single sterile barrier system with protective packaging inside
Glucose	Use glucose in the purge fluid
 Do Not Flush	Do Not Flush with syringe
Do Not Clean with Alcohol	Do not use alcohol or alcohol-based products for cleaning.

IMPELLA 5.5 CATHETER PARAMETERS

Speed range	0 to 33,000 rpm
Power consumption	Less than 16 W
Voltage	Max. 20 V DC
Flow-Maximum	5.5 L/min
Purging the Impella 5.5 Catheter	5% glucose solution in water with heparin concentration of 25 or 50 IU per mL
Recommended purge fluid	
Glucose concentration	5% to 20%
Purge pressure	300 to 1100 mmHg
Purge flow	2 to 28 mL/hr
Maximum duration of use	
EU	Up to 30 days
Reliability	91.4% with 80% confidence at 30 days based on chi-squared test
Dimensions of Impella 5.5 Catheter	
Length of invasive portion (without catheter)	114 mm
Diameter (Motor)	18 Fr (6.3mm)
Diameter (Cannula)	21 Fr (7 mm)
Classification per DIN EN 60601-1	Protection class I, degree of protection: CF (Automated Impella Controller and Impella 5.5 Catheter)
Classification per directive 93/42/EEC	Class III
Latex free	Yes

IMPELLA 5.5® WITH SMARTASSIST® CATHETER DIMENSIONS

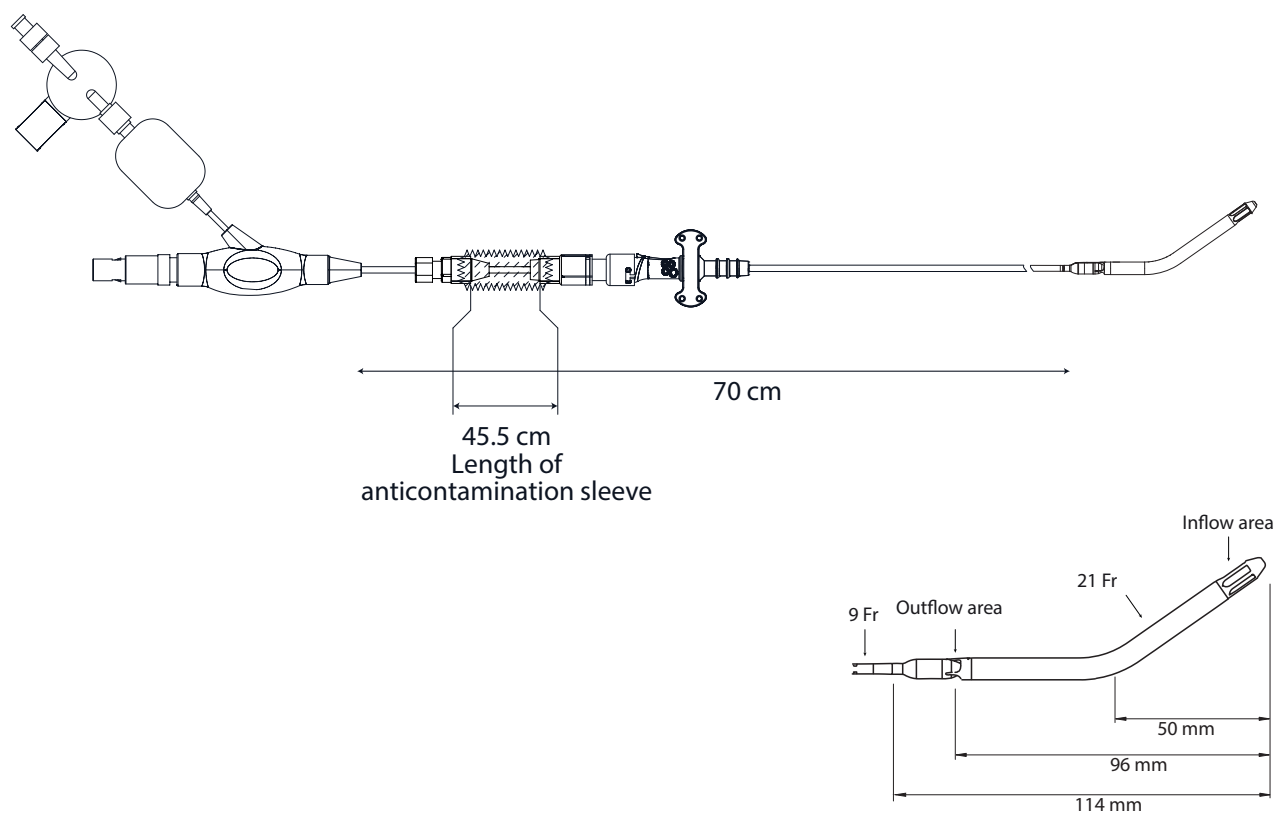DISPOSING OF THE IMPELLA 5.5 WITH SMARTASSIST CATHETER AND ACCESSORIES (EUROPEAN UNION)

The Impella 5.5 with SmartAssist Catheter is disposable and must be disposed of in accordance with hospital regulations for blood contaminated materials.

The Automated Impella Controller is marked according to Directive 2012/19/EU. Devices sold within the EU can be returned to Abiomed Europe GmbH for correct disposal.

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Abiomed, Inc.